

Floorstanding Cryostat MNT

The fully automatic floorstanding model
with optional extra features

SLEE
M A I N Z

SLEE
M A I N Z

Weight and Dimensions	Depth x Width x Height	
Dimensions (in mm)	760 x 750 x 1140	
Dimensions packed (in mm)	880 x 840 x 1330	
Weight net	110 kg	
Weight packed	170 kg	
Electrical parameters	Art.-No. 1015000	230 V / 50 Hz
	Art.-No. 1015100	230 V / 60 Hz
	Art.-No. 1015200	115 V / 60 Hz
	Art.-No. 1015300	100 V / 50 Hz
	Art.-No. 1015400	100 V / 60 Hz

Please ask for other units for Histopathology and consumables as well as laminar flow equipment.



Benchtop Cryostat MTC –
optional with height adjustable support



Floorstanding Cryostat MEV

(Specifications are subject to change without notice)

Floorstanding Cryostat MNT

The fully automatic floorstanding model
with optional extra features



Microtomes

Cryostats

MTC

MEV

► MNT

Paraffin

Embedding Centers

Tissue processors

Stainers

Cover slipper

Knife sharpener

Floating bath

Slide warmer

SLEE MEDICAL GMBH
P. O. Box 10 05 07
D-55129 Mainz · Germany
Fon +49 (0) 6131/95871-0
Fax +49 (0) 6131/95871-0
mail@sleemainz.com
www.sleemainz.com

Floorstanding Cryostat MNT

The fully automatic floorstanding model with optional extra features

The model MNT has a rotary microtome with motorised advance / return and cutting stroke.

- ▶ Modern rotary microtome with motorised μm advance/return and coarse advance for use of standard knife holder or disposable blade holder (exchangeable)
- ▶ Cutting motor with adjustable speed
- ▶ New anti-roll unit with large anti-roll plate
- ▶ Viewing window consists of dual insulation glass, frost-free
- ▶ Clear electronic control for chamber temperature, automatic defrost and LED-light, with real time clock, with power-saving weekend mode
- ▶ Fully adjustable automatic defrost, several times per day possible
- ▶ Low temperature storage area for 24 specimen holders
- ▶ Optional:
 - Shock-freezing to -60°C
 - Peltier shock freezing to -55°C (2 positions)
 - Specimen cooling to -45°C
 - Debris extraction
 - Automatic disinfection system on thermal and chemical basis
 - UV-C disinfection


SLEE
M A I N Z



1. Electronic control with foil pad operation and LCD display
2. Key control for section thickness adjustment etc.
3. Compact, reliable cooling system with CFC-free cooling gas R 404 A
4. Handwheel
5. Frost-free insulation glass
6. Rim-heating system for insulation glass
7. Hand wheel lock, ideally positioned
8. Large chamber with drain for hygienic cleaning

Technical Data

Cooling capacity	to -40°C
Section thickness range	0 – 2 μm in 0,5 μm steps 2 – 20 μm in 1 μm steps 20 – 50 μm in 2 μm steps 50 – 100 μm in 5 μm steps
Trimming adjustable	from 0 – 300 μm , free programmable
Specimen holder	22 or 35 mm dia. with rubber surround for handling comfort
Specimen support	orientable on 3 planes, XY-Axes 8° , Z-Axes 360°
Vertical stroke	57 mm
Horizontal specimen advance	28 mm

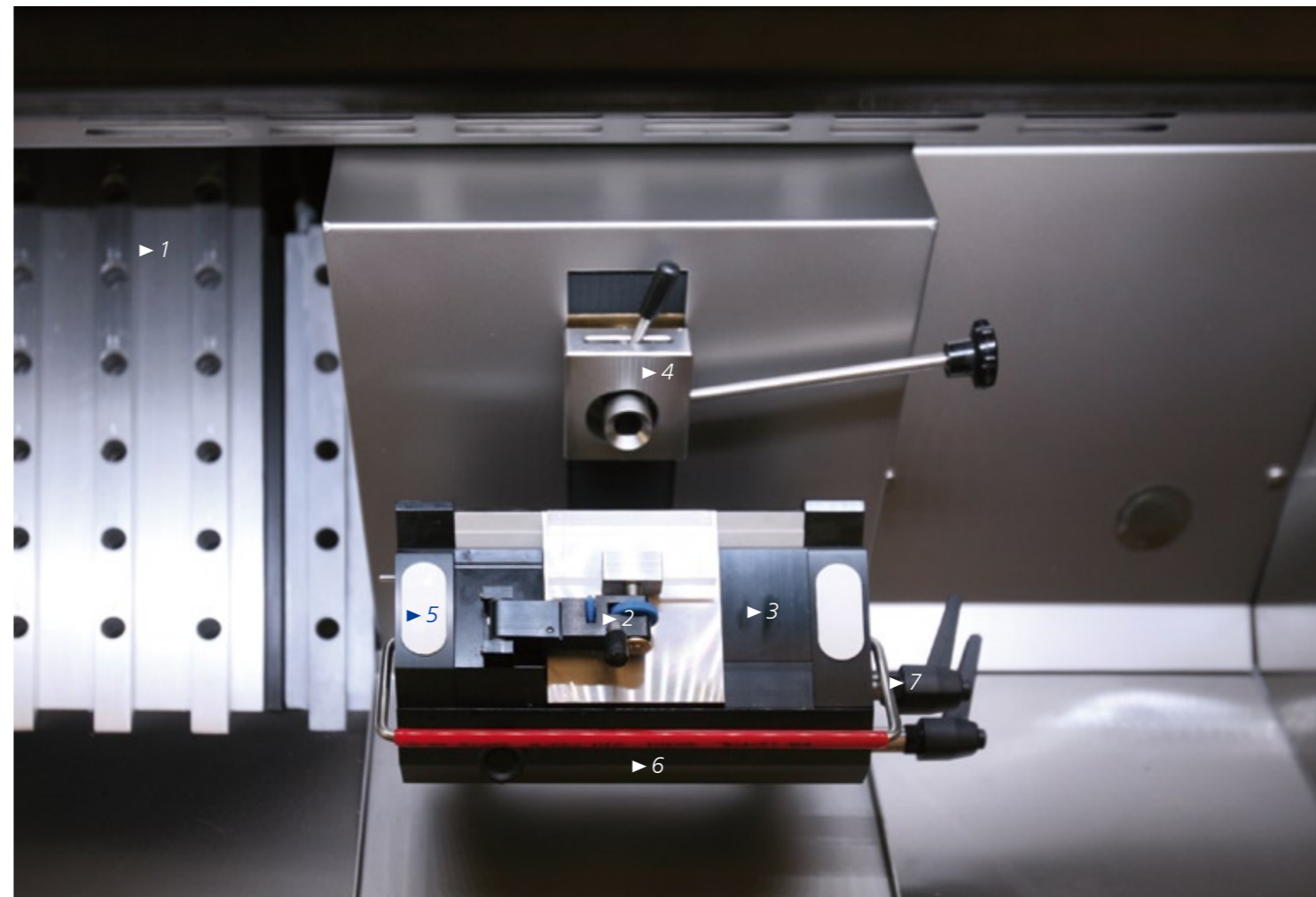

1. 24 freezing positions to -45°C
2. Anti-roll assembly
3. Standard knife holder (for steel- or hard metal knives) or disposable blade holder
4. Object orientation: XY Axes 8° , Z-Axes 360°
Specimen holder 22 or 35 mm dia. with rubber surround for handling comfort
5. Plastic freeze protection for comfortable operation
6. Knife holder base
7. Integrated finger protection