


Cover slipper MCS

One of the worldwide fastest
automated units for efficient cover slipping

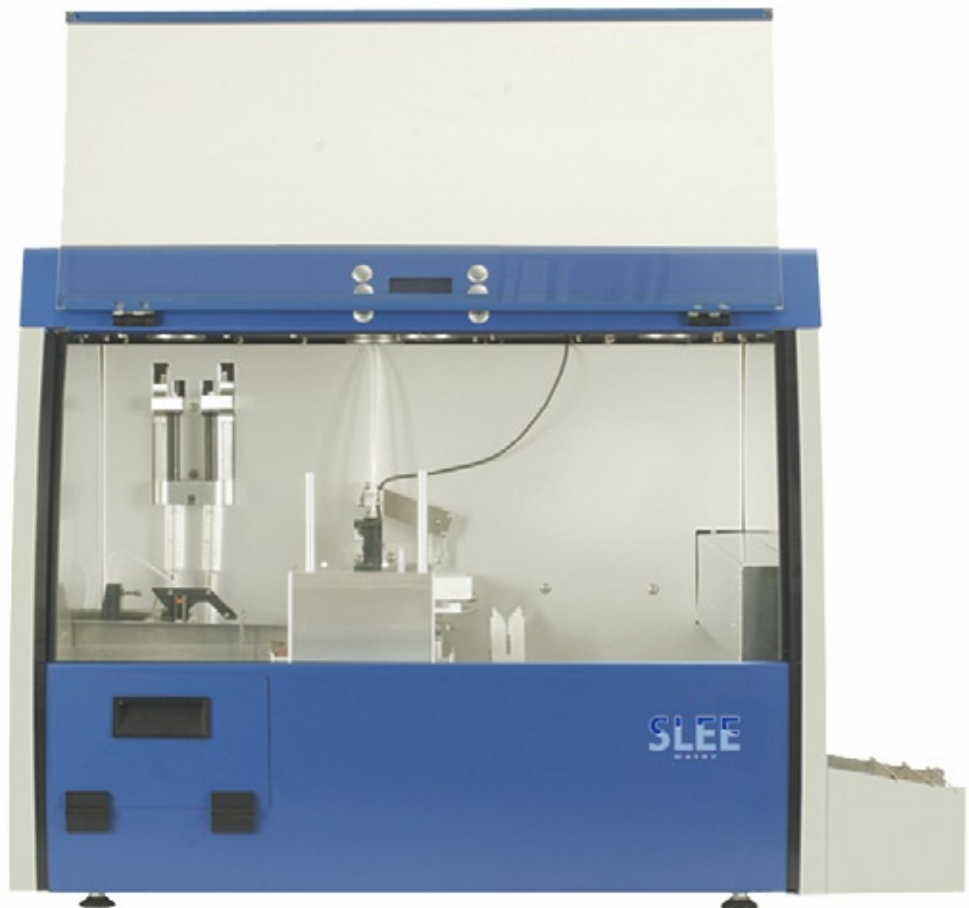

Microtomes

Cryostats

Paraffin

Embedding Centers

Tissue processors

Stainers

Cover Slipper

► **MCS**

Knife sharpener

Floating bath

Slide warmer

Cover slipper MCS

One of the worldwide fastest automated units for efficient cover slipping



The innovative cover slipping technology allows fast identification of single stainings as well as clear classification. The MCS works with highest quality standards, cover slipping is done for each single glass separately fast and without bubbles. This method guarantees a reliable throughput.

The cover slipper MCS can be ideally used in combination with the stainer MAS. The modular design simplifies the step of manual hand over of the sample holders to full automation – also in times of small budgets. The combination uses a minimum of space.

- ▶ Capacity 120 pieces continuously up to 520 slides per hour
- ▶ Suitable for commercial cover slipping medium containers
- ▶ Universal power input 100 – 240 V, 50/60 Hz
- ▶ Active carbon filter for user safety



Stainer MAS in combination with Cover Slipper MCS



Stainer MSM for universal use

Technical Data

Cover slipping time	7 – 12 sec. per slide
Cover slipping amount	40 – 250 µl, adjustable

Specimen glasses

Dimensions standard	24 x 40, 50, 55, 60 mm
Dimensions optional	24 x 32, 36 mm, 25 x 40, 50, 55, 60 mm
Thickness	0,12 – 0,17 mm
Capacity	350 pcs.

Weight and Dimensions

Dimensions (in mm)	670 x 900 x 640
Dimensions packed (in mm)	800 x 1200 x 800

Weight net	75 kg
Weight packed	85 kg

Electrical parameters	Art.-No. 14000000	230 V / 50 Hz
	Art.-No. 14000001	230 V / 60 Hz
	Art.-No. 14000002	115 V / 60 Hz
	Art.-No. 14000003	100 V / 50 Hz
	Art.-No. 14000004	100 V / 60 Hz

Please ask for other units for Histopathology and consumables as well as laminar flow equipment

(Specifications are subject to change without notice)

SLEE MEDICAL GMBH

P. O. Box 10 05 07
D-55129 Mainz · Germany

Fon +49 (0) 6131/95871-0
Fax +49 (0) 6131/95871-0

mail@sleemainz.com
www.sleemainz.com